

Prepared for:

**Centers for Medicare & Medicaid Services and Office of the
National Coordinator for Health Information Technology**

**CMS Alliance to Modernize Healthcare
Federally Funded Research and Development Center**

**Technical Authority for the Unified Clinical Quality
Improvement Framework**

Bonnie User Guide

Version 3.2

May 21, 2019

The views, opinions, and/or findings contained in this report are those of The MITRE Corporation and should not be construed as official government position, policy, or decision unless so designated by other documentation.

© 2019, The MITRE Corporation. All Rights Reserved. Approved for Public Release. Distribution Unlimited 13-4297.

Record of Changes

Version	Date	Author / Owner	Description of Change
0.1	April 15, 2014	Andre Quina / MITRE	Initial draft
0.1.1	July 18, 2014	Cristen Jones / MITRE	Update for Bonnie v 1.1
0.2	April 13, 2015	Cristen Jones / MITRE	Update for Bonnie v 1.3
1.4	January 13, 2016	Pace Ricciardelli / MITRE	Update for Bonnie v 1.4
1.5	November 18, 2016	David Connolly and Lizzie DeYoung / MITRE	Update for Bonnie v 1.5
1.6	January 4, 2017	Kristian Mulcahy and Chris Tohline / MITRE	Update for Bonnie v 1.6
1.6.5	June 8, 2017	Chris Tohline / MITRE	Update for Bonnie v 1.6.5
2.0	November 14, 2017	Rose Almonte, Lizzie Charbonneau, Chris Tohline, Jenn Reed / MITRE	Update in preparation for Bonnie v 2.0
2.1	May 3, 2018	Lizzie Charbonneau / MITRE	Update for Bonnie v 2.1
2.1.2	May 31, 2018	Lizzie Charbonneau / MITRE	Update for Bonnie v 2.1.2
2.2	July 13, 2018	Cole Springate-Combs / MITRE	Update for Bonnie v 2.2
2.3	August 23, 2018	James Bradley / MITRE	Update for Bonnie v 2.3
3.0	October 16, 2018	James Bradley / MITRE	Update for Bonnie v 3.0
3.1	November 15, 2018	James Bradley / MITRE	Update for Bonnie v 3.1
3.2	May 21, 2019	James Bradley / MITRE	Update for Bonnie v 3.2

Table of Contents

1. Introduction	1
1.1 Background	1
1.2 Purpose	1
1.3 Application Description	1
2. User Account Creation	3
2.1 Login Page.....	3
2.2 Creating a New User	3
2.3 Resetting a Password.....	4
2.4 Account Management	5
3. Measure Dashboard.....	6
3.1 Overview	6
3.2 Loading a New Measure	7
3.2.1 Value Set Options	9
3.3 Updating an eCQM	10
3.4 Creating Synthetic Test Records	11
4. Measure Results View.....	12
4.1 Overview	12
4.2 Measure Logic.....	14
4.3 Creating a New Test Record	14
4.4 Calculation Results.....	15
4.5 Editing a Test Record.....	17
4.6 Cloning a Test Record.....	17
4.7 Deleting a Test Record.....	17
4.8 Updating a Measure	18
4.9 Deleting a Measure.....	18
5. Building a Patient Test Record.....	19
5.1 Overview	19
5.2 Building a Synthetic Patient.....	20
5.3 Building the Patient History	21
5.3.1 Patient History Items that Are Related to Past Items.....	23
5.4 Incremental Calculation	23
6. Feedback and Support.....	24
7. Frequently Asked Questions	25
Acronyms.....	27

List of Figures

Figure 1. Bonnie Login Page	3
Figure 2. Account Registration Page	4
Figure 3. Password Reset Page	4
Figure 4. Account Management Page	5
Figure 5. Measure Dashboard View	6
Figure 6. New Measure Dialog	7
Figure 7. Finalize Measure Dialog	8
Figure 8. Value Set Options – Profile	9
Figure 9. Value Set Options - Release	10
Figure 10. Updating Measure Dialog	11
Figure 11. Measure View	13
Figure 12. Measure Terminology and Overlapping Value Sets	14
Figure 13. Expanded Results View	15
Figure 14. Logic Calculation Highlight – Passing Results	16
Figure 15. Logic Calculation Highlight – Expanded Result	16
Figure 16. Logic Calculation Highlight – Failing Results	17
Figure 17. Patient Builder View	19
Figure 18. Continuous Variable Measures Expected Populations	21
Figure 19. Building Patient History, Including Edit Clinical Element View	22
Figure 20. References Section of the Patient History Builder	23
Figure 21. User Group Link on Bonnie Splash Page	24
Figure 22. User Group Link in the Application Header	24

1. Introduction

1.1 Background

Bonnie is a software tool that allows electronic clinical quality measure (eCQM) developers to test and verify the behavior of their eCQM logic. The main goal of the Bonnie application is to reduce the number of defects in eCQMs by providing a robust and automated testing framework. The Bonnie application allows measure developers to independently load measures that they have constructed using the Measure Authoring Tool (MAT). Developers can then use the measure metadata to build a synthetic patient test deck for the measure from the clinical elements defined during the measure construction process. By using measure metadata as a basis for building synthetic patients, developers can quickly and efficiently create a test deck for a measure. The Bonnie application helps measure developers execute the measure logic against the constructed patient test deck and evaluate whether the logic aligns with the intent of the measure.

Bonnie has been designed to integrate with the nationally recognized data standards the Centers for Medicare & Medicaid Services (CMS) quality reporting programs use for expressing eCQM logic for machine-to-machine interoperability. This integration provides enormous value to the eCQM program and federal policy leaders and stakeholders. The Bonnie tool verifies that the new and evolving standards for eCQMs used in the CMS quality reporting programs are flexible and can be implemented in software.

Bonnie was also designed to provide an intuitive and easy-to-use interface based on feedback from the broader measure developer community. A key goal of the Bonnie application is to deliver a user experience that provides an efficient and intuitive method for constructing synthetic patient records for testing and validating eCQMs.

Finally, the Bonnie software is freely available via an Apache 2.0 open source license. The CMS Electronic Health Records (EHR) Incentive Program makes all or parts of the Bonnie software available for inspection, verification, and even reuse by other government programs or federal contractors.

1.2 Purpose

The purpose of this document is to describe the functionality of the Bonnie web application that allows measure developers to test and verify the behavior of their eCQM logic. This document provides Bonnie users with step-by-step instructions for testing eCQMs by building synthetic patient records.

1.3 Application Description

The Bonnie application provides the capability to import measures defined in Health Quality Measure Format (HQMF) Extensible Markup Language (XML) and Clinical Quality Language (CQL). The HQMF specification provides the metadata and the CQL libraries provide the logic for calculating an eCQM. The Bonnie application can load an eCQM containing the HQMF and CQL describing a measure and programmatically convert both into an executable format that allows calculation of the measure directly from the specification.

The CMS MAT is the primary source for eCQMs used by the Bonnie application. Measure developers use the MAT to build eCQMs and export those measures as measure bundles containing HQMF and one or more CQL libraries. These measure bundles can be downloaded from the MAT and loaded into the Bonnie user interface for measure testing.

Once an eCQM has been loaded into the Bonnie application, a user can inspect the measure logic and then build synthetic test records and set expectations on how those test records will calculate against a measure. This capability to build synthetic test patient records, set expectations against those records, and calculate the measures using those patient records provides an automated and efficient testing framework for eCQMs. Through the Bonnie-supported eCQM testing framework, measure developers can more clearly understand the behavior of the measure logic and validate that the measure logic encodes their intent, as well as validate multiple iterations of measure updates against a test deck.

In addition, the development of a test deck as part of measure development provides benefits after the measures are finalized. For example, the test deck can demonstrate the intent of the measure through the patient examples in the deck. The systems in the test deck implement those measures with a means to validate their development through a base set of synthetic patient records with known expectations for calculating against the implemented measures. Finally, the test deck can provide input for the Office of the National Coordinator for Health Information Technology (ONC) Certification Program. If the synthetic patients are clinically valid, Cypress¹ users can use them in addition to or instead of their own patients.

¹ Cypress is the rigorous and repeatable testing tool of Electronic Health Records (EHR) and EHR modules used in calculating electronic clinical quality measures used in CMS's quality reporting programs.

2. User Account Creation

2.1 Login Page

Users require a valid account to access the Bonnie application. All eCQMs loaded into the Bonnie application and all synthetic test patients are isolated by account. Therefore, users can only view, access, and modify data that they have loaded under their own account. Figure 1 shows the login screen for the Bonnie application. To log in, a user must provide the email address and password for a valid account.


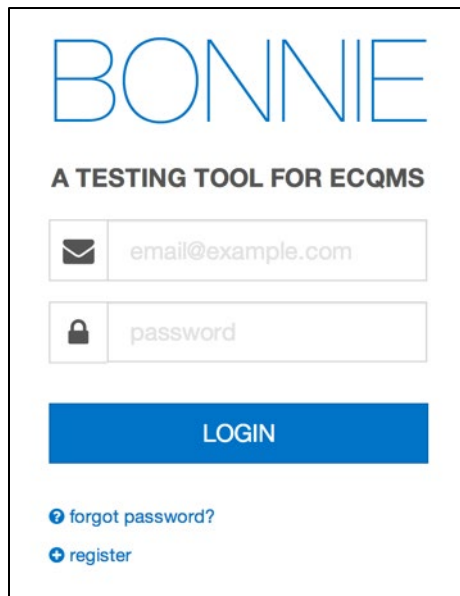
The image shows the Bonnie login page. At the top, the word "BONNIE" is displayed in a large, blue, sans-serif font. Below it, the text "A TESTING TOOL FOR ECQMS" is written in a smaller, bold, black font. There are two input fields: the first is for an email address, with a placeholder "email@example.com" and an envelope icon to its left; the second is for a password, with a placeholder "password" and a lock icon to its left. Below these fields is a blue rectangular button with the word "LOGIN" in white, uppercase letters. At the bottom of the form, there are two links: "forgot password?" with a question mark icon and "register" with a plus icon.

Figure 1. Bonnie Login Page

2.2 Creating a New User

A user can create a new account by clicking the “register” link on the login page. The register link brings the user to account creation page shown in Figure 2. A user can create a new account by filling out the fields in the registration form and clicking the “REGISTER” button. Once an account has been created, the user can log in to the Bonnie application with the email address and password specified as part of account creation.

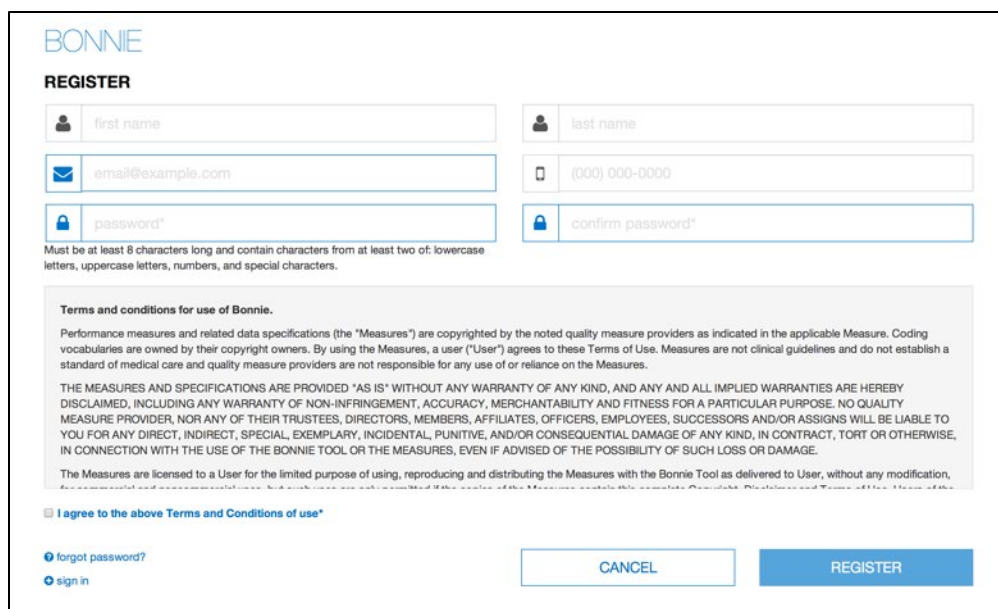

BONNIE

REGISTER

first name

last name

email@example.com

(000) 000-0000

password*

confirm password*

Must be at least 8 characters long and contain characters from at least two of: lowercase letters, uppercase letters, numbers, and special characters.

Terms and conditions for use of Bonnie.

Performance measures and related data specifications (the "Measures") are copyrighted by the noted quality measure providers as indicated in the applicable Measure. Coding vocabularies are owned by their copyright owners. By using the Measures, a user ("User") agrees to these Terms of Use. Measures are not clinical guidelines and do not establish a standard of medical care and quality measure providers are not responsible for any use of or reliance on the Measures.

THE MEASURES AND SPECIFICATIONS ARE PROVIDED 'AS IS' WITHOUT ANY WARRANTY OF ANY KIND, AND ANY AND ALL IMPLIED WARRANTIES ARE HEREBY DISCLAIMED, INCLUDING ANY WARRANTY OF NON-INFRINGEMENT, ACCURACY, MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. NO QUALITY MEASURE PROVIDER, NOR ANY OF THEIR TRUSTEES, DIRECTORS, MEMBERS, AFFILIATES, OFFICERS, EMPLOYEES, SUCCESSORS AND/OR ASSIGNS WILL BE LIABLE TO YOU FOR ANY DIRECT, INDIRECT, SPECIAL, EXEMPLARY, INCIDENTAL, PUNITIVE, AND/OR CONSEQUENTIAL DAMAGE OF ANY KIND, IN CONTRACT, TORT OR OTHERWISE, IN CONNECTION WITH THE USE OF THE BONNIE TOOL OR THE MEASURES, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH LOSS OR DAMAGE.

The Measures are licensed to a User for the limited purpose of using, reproducing and distributing the Measures with the Bonnie Tool as delivered to User, without any modification, for non-commercial use. Such use is not permitted if the use of the Measures results in the creation of a new product or service.

☐ I agree to the above Terms and Conditions of use*

[forgot password?](#)


[sign in](#)

CANCEL REGISTER

Figure 2. Account Registration Page

2.3 Resetting a Password

If a user forgets a password or an account is locked, the user can reset the password using the password reset page shown in Figure 3. This page is accessed from the “forgot password?” link on the login page (Figure 1). On the password reset page, the user provides the email address associated with the account and then presses the “Send” button. This action sends an email to the registered email address for the account to facilitate resetting the user’s account password.

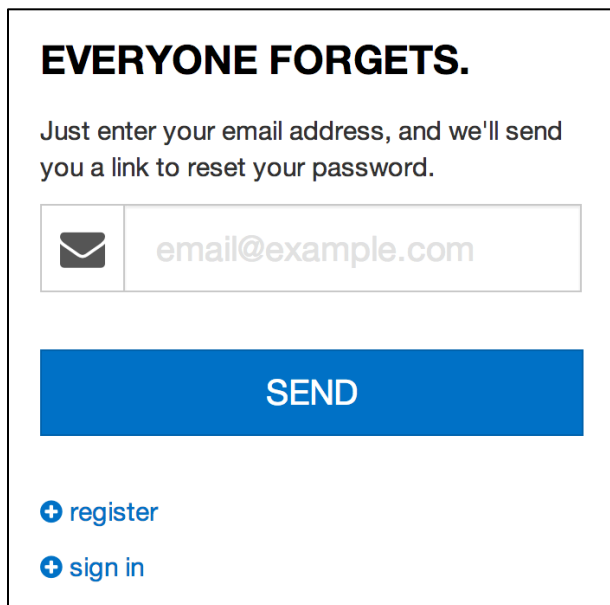

EVERYONE FORGETS.

Just enter your email address, and we'll send you a link to reset your password.

email@example.com

SEND


[+ register](#)

[+ sign in](#)

Figure 3. Password Reset Page

2.4 Account Management

After logging into the application, the user can change the information associated with the account by accessing the account management page shown in Figure 4. The account management page can be opened by clicking the “Account” link in the application header shown at the top of Figure 4. From this page, users can change the information provided during the registration process and select a new password for their account.

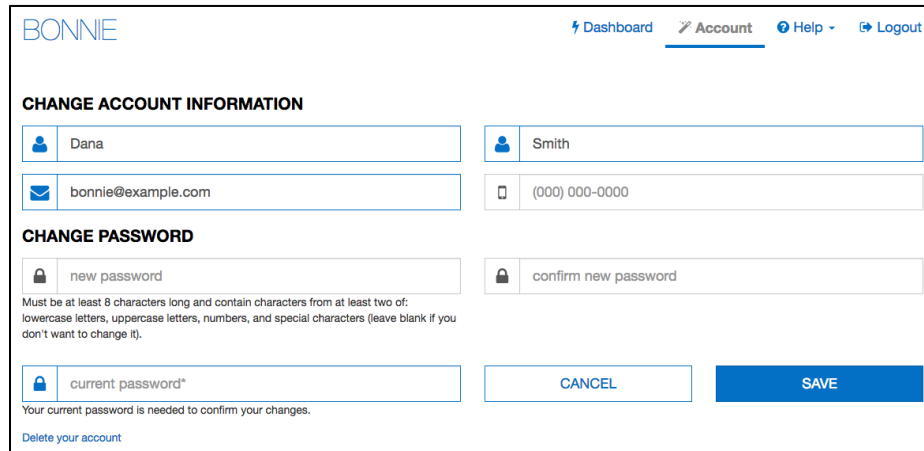

The screenshot displays the 'BONNIE' application interface. At the top, a navigation bar includes links for 'Dashboard', 'Account' (which is the active page), 'Help', and 'Logout'. The main content area is titled 'CHANGE ACCOUNT INFORMATION' and contains two columns of input fields. The left column includes fields for 'First Name' (containing 'Dana'), 'Email' (containing 'bonnie@example.com'), 'New Password' (containing 'new password'), and 'Current Password' (containing 'current password*'). The right column includes fields for 'Last Name' (containing 'Smith') and 'Phone Number' (containing '(000) 000-0000'). Below the 'New Password' field, a note specifies that the password must be at least 8 characters long and contain characters from at least two of the following categories: lowercase letters, uppercase letters, numbers, and special characters. At the bottom of the form, there are 'CANCEL' and 'SAVE' buttons. A link labeled 'Delete your account' is located at the bottom left of the page.

Figure 4. Account Management Page

3. Measure Dashboard

3.1 Overview

The Measure Dashboard page, as shown in Figure 5, is the initial page presented to users when they log into the application. The Measure Dashboard shows the calculation status of each eCQM loaded into the system. The calculation status shows how many patient records have been built for the measure, whether the measure is currently passing or failing, and how many patients are passing or failing for each measure.

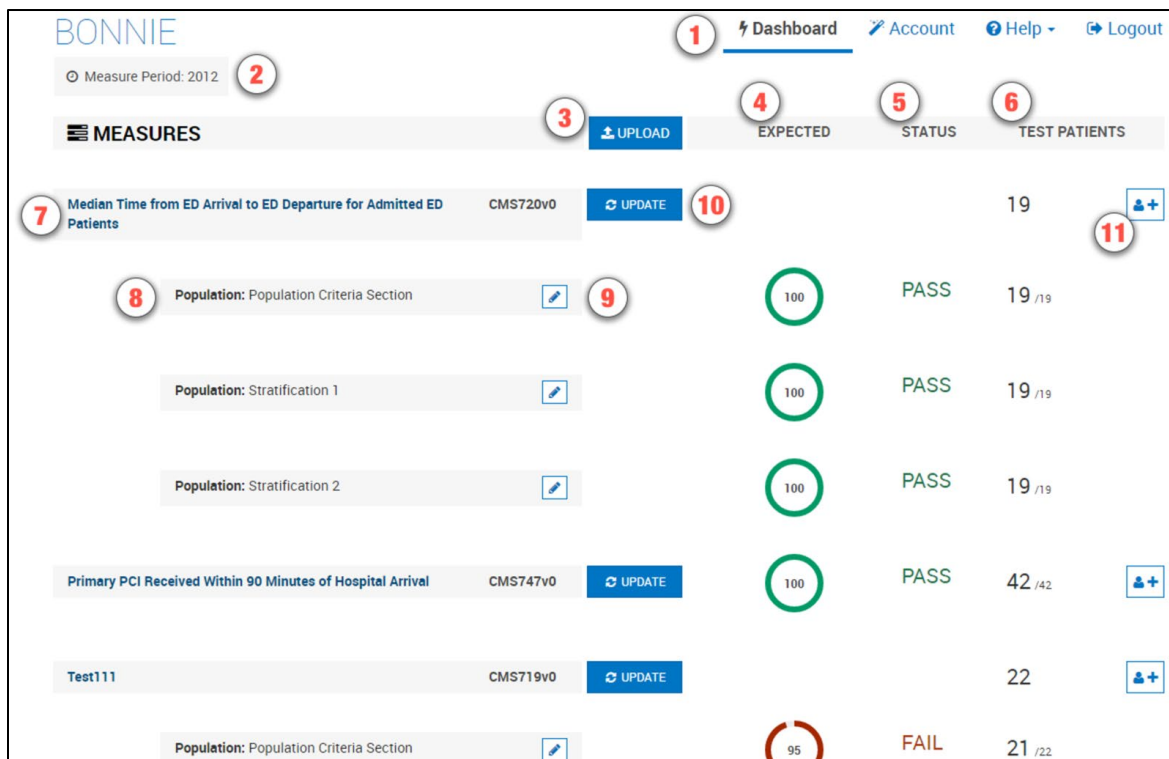


Figure 5. Measure Dashboard View

The Measure Dashboard View employs the following user interface (UI) elements (indicated by their item numbers in Figure 5):

1. Header – Allows the user to navigate to different parts of the application, access account information, send a support email (Contact), and log out of the application.
2. Measure Period Date – Displays the measurement period used for calculating eCQMs.
3. Upload Button – Allows the user to upload a new eCQM.
4. Expected Column – Displays the percentage of passing patients for the eCQM.
5. Status Column – Displays the status of the eCQM (New, Pass, Fail).
6. Test Patient Column – Displays the number of patients passing out of the total number of patients.

7. Measure Title – Displays the title for the eCQM and allows navigation to the Measure Dashboard View.
8. Subpopulation and stratification titles – Displays the titles for subpopulations or stratifications of an eCQM.
9. Edit Title Button – Allows the user to rename a subpopulation or stratification.
10. Update Button – Allows the user to update a previously loaded eCQM.
11. Add Patient Button – Allows the user to start building a new patient for an eCQM.

3.2 Loading a New Measure

When users log into the system for the first time, no eCQMs will be associated with the account. The user's first step is to load an eCQM into the account to begin testing the eCQM with the Bonnie application. The New Measure dialog, as shown in Figure 6, prompts the user to upload an eCQM.

The user may either upload an eCQM now or close this display and upload an eCQM later. After one or more eCQMs have been uploaded, this page displays the current set of eCQMs loaded into the system along with the subpopulations and stratifications associated with the eCQMs. The Measure Dashboard allows users to navigate to the details of individual eCQMs, upload a new eCQM, or update the definition of an existing eCQM. To upload a new eCQM, click the “Upload” button (item #3 in Figure 5) on the Measure Dashboard, which opens the New Measure Dialog shown in Figure 6.

The screenshot shows the 'New Measure' dialog box with the following elements and numbered callouts:

- 1** MEASURE DATA: A blue 'BROWSE' button.
- 2** VSAC USERNAME: A text input field containing 'username'.
- 3** VSAC PASSWORD: A text input field containing 'password'.
- 4** Register for VSAC: A blue link with a plus icon.
- 5** VALUE SETS: Two buttons, 'Profile' (selected) and 'Release'.
- 6** Latest eCQM: A dropdown menu showing 'Latest eCQM (eCQM Update 2018-09-17)'.
- 7** Draft: A checkbox that is checked.
- 8** Override options with measure defined versions: A checkbox that is unchecked.
- 9** TYPE: Two buttons, 'Eligible Professional' (selected) and 'Eligible Hospital'.
- 10** CALCULATION: Two buttons, 'Patient Based' (selected) and 'Episode of Care'.
- 11** Include Supplemental Data Element Calculations: A checkbox that is unchecked.
- 12** CLOSE: A blue button.
- 13** LOAD: A blue button.

Figure 6. New Measure Dialog

After clicking the “Upload” button, the user needs to fill out the information requested from the New Measure Dialog. These fields are outlined below and are indicated by their corresponding item numbers in Figure 6.

1. Browse – Opens a file browser enabling the user to select a MAT measure package to upload into Bonnie.
2. VSAC Username and Password – Fields for entering the user’s VSAC credentials to enable Bonnie to load the value sets included within the eCQM.
3. Value Set Options – Options to define which value sets should be downloaded from VSAC. More information about these options is included in subsection 3.2.1.
4. Type – Option to define if the eCQM is an Eligible Professional/Eligible Clinician or Eligible Hospital/Critical Access Hospital eCQM.
5. Calculation – Option to define if the eCQM is a Patient-based eCQM or an Episode of Care-based eCQM.
6. Include Supplemental Data Element Calculations – Option to define if supplemental data elements should be included in calculations and displayed in the supplemental data element section in the eCQM logic.
7. Close – Closes the New Measure Dialog without uploading the eCQM.
8. Load – Uploads the specified eCQM using the selected options.

The action of clicking the “Load” button in the New Measure Dialog uploads the eCQM to the application for processing. If the eCQM loaded is Episode of Care or has multiple populations, the user is presented with the Finalize Measure Dialog shown in Figure 7. This dialog allows the user to specify the episode(s) of care for the eCQM and provide titles for subpopulations. After filling out the Finalize Measure fields, the user clicks the “Done” button to finish loading the eCQM. Once eCQM loading is complete, the application directs the user to the Measure Dashboard (Figure 5) with the new measure available. If the measure is not Episode of Care or does not have subpopulations, the Finalize Measure Dialog is not presented because no additional information is required to load the measure.

Figure 7. Finalize Measure Dialog

3.2.1 Value Set Options

Bonnie allows users to download value sets associated with an eCQM based on either a VSAC profile or VSAC release. Bonnie also allows users to download draft value sets (if “Profile” is selected) and/or override their selections with value set versions specified within the CQL logic itself as shown in Figure 8.

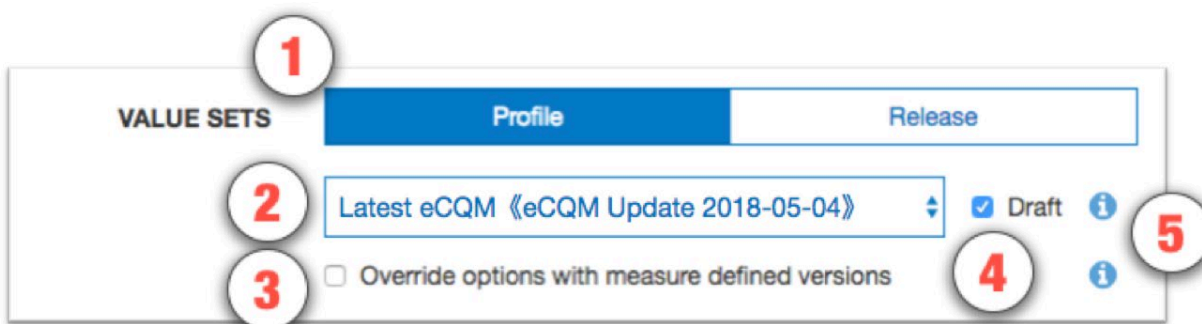


Figure 8. Value Set Options – Profile

Figure 8 shows the Value Set Options section (item #3 in Figure 6) with the default selections. The default selections are “Profile” with the “Latest eCQM” profile selected and “Draft” checked. These default selections have the same behavior as the default selections in Bonnie v2.0. The UI elements included when “Profile” is selected are as follows:

1. **Profile** – When Profile is selected, Bonnie queries the Value Set Authority Center (VSAC) for value sets calculated by an expansion profile. Expansion profiles are instructions for calculating value set content based on specific code system versions and retired legacy codes.
2. **Profile dropdown selector** – Allows the user to select an expansion profile to use for retrieving value sets. The default selection is “Latest eCQM <<{profile}>>” where {profile} is the VSAC profile referenced by “Latest eCQM”. In this screen shot, the referenced profile is “eCQM Update 2018-05-04”, so the display shows “Latest eCQM <<eCQM Update 2018-05-04>>”.
3. **Override options with eCQM-defined versions** – It is possible for an eCQM to contain value set versions within the CQL itself. This option allows the user to override the selected option to use the eCQM-defined value set versions. If this option is selected and a version for a value set is **not** defined within the measure, then Bonnie will use the selected profile to retrieve the appropriate value set version.
4. **Draft** – When this option is selected, Bonnie will retrieve the draft value set expansions for the selected expansion profile. Note that if the user selects draft and is not an assigned author or steward group member for one or more of the value sets included in the measure, the measure upload will fail. Draft is selected by default.
5. **Information buttons** – Contains additional information regarding the value set options.

Figure 9 shows the Value Set Options section (item #3 in Figure 6) with “Release” selected.

The screenshot shows the 'VALUE SETS' section with the 'Release' tab selected. The 'Profile' tab is also visible. Below the tabs, there is a dropdown menu labeled 'CMS eCQM' (2) and another dropdown menu labeled 'eCQM Update 2018-05-04' (3). Below these, there is a checkbox labeled 'Override options with measure defined versions' (4). To the right of the dropdowns, there are information buttons (5).

Figure 9. Value Set Options - Release

The UI elements included when “Release” is selected are as follows:

1. Release – When Release is selected, Bonnie queries VSAC for the published value sets associated with the specified release.
2. Program dropdown selector – Allows the user to select the program that includes the release they would like to use.
3. Release dropdown selector – Allows the user to select the appropriate release.
4. Override options with measure-defined versions – It is possible for an eCQM to contain value set versions within the CQL itself. This option allows the user to override the selected option to use the measure-defined value set versions. If this option is selected and a version for a value set is **not** defined within the eCQM, then Bonnie will use the selected release to retrieve the appropriate value set version.
5. Information buttons – Contains additional information regarding the value set options.

3.3 Updating an eCQM

Once the eCQM has been loaded, the testing process may identify issues with the eCQM. When issues are identified, the logic must be updated in the MAT to resolve these issues. Alternatively, the measure could be updated in the MAT as part of an Annual Update. After an eCQM has been updated in the MAT, it may be necessary to update that eCQM in Bonnie for testing. To update an eCQM, follow these steps:

1. Click the “Update” button (item #10 in Figure 5) on the Measure Dashboard, which displays the “Update Measure” dialog as shown in Figure 10.
2. Select a new MAT export zip with the updated eCQM definition.
3. Click the “Load” button to load the new version of the eCQM.

Test104

MEASURE DATA **BROWSE** DRAFT CMS104 CQL_with_json.zip

VSAC USERNAME

VSAC PASSWORD
[Register for VSAC](#)

VALUE SETS **Draft** Measure Defined

TYPE Eligible Professional (EP)

CALCULATION Patient Based

REBUILD PATIENTS **No** Yes

CLOSE **LOAD**

Figure 10. Updating Measure Dialog

3.4 Creating Synthetic Test Records

Once a set of eCQMs has been loaded into the Bonnie application, users can start building test patients for the eCQMs. To build a test patient from the Measure Dashboard (Figure 5), click the “Add Patient” button (item #12). This action opens the patient builder screen as shown in Figure 17. For more information, please refer to Section 5, *Building a Patient Test Record*.

After the user creates synthetic test patients for eCQMs, the Measure Dashboard displays summary calculation results of the patients associated with each eCQM loaded by the user (Figure 5). As shown in Figure 5, UI elements numbered 5, 6, and 7 on the Measure Dashboard provide the summary results for each eCQM. In the “Expected” column (item #5), the UI displays the percentage of patients associated with the eCQM whose calculated values meet the expectations set for the patient. The “Status” column (item #6) indicates the current state of the measure—whether it is passing, failing, or new. An eCQM is in the passing state if all associated patients are meeting expectations. If one or more patients are not meeting expectations, then the eCQM is in a failing state. Any eCQM is in the new state when there are no patients associated with it. Finally, the “Test Patients” column (#7) displays the number of passing patients over the total number of patients as a fraction. The measure view, as represented in Figure 11 in Section 4, shows more detailed results for an eCQM.

4. Measure Results View

4.1 Overview

As shown in Figure 11, the Measure View page displays the detailed information, associated patients, and calculation results for a single eCQM. In this view, users can add new patients to an eCQM, update an eCQM, and delete an eCQM. To access the Measure View, click on the “Measure Title” link (item #7) on the Measure Dashboard, as shown in Figure 5.

The Measure View page presents the title and description of the eCQM, the logic for the eCQM, and a measure complexity indicator. The page lists the current set of test patients associated with the eCQM in the column on the right-hand side of the page. This column displays the summary calculation results for the patient test deck associated with the eCQM along with the list of patients and the individual calculation results for each. This page also shows the results of calculation for a single patient along with an overlay of the calculation results on the eCQM logic.

The following UI elements appear at the top of the Measure View page (indicated by their item numbers in Figure 11):

1. Measure Title – Displays the title and description of the eCQM.
2. Measure Actions – Allows the user to delete or update an eCQM definition.
3. Measure Complexity – Indicates the eCQM’s degree of complexity.
4. Measure Subpopulations or Stratifications – Allows access to different subpopulations or stratifications in the eCQM.
5. Measure Logic – Displays a representation of the logic for the eCQM.
6. Logic Highlighted With Coverage – Displays the logic for the eCQM, highlighting which lines of the eCQM are covered by the test patients.
7. Patient Actions – Allows the user to export patient records as QRDA or Excel.
8. Percent Successful – Displays the percent of patients currently meeting expectations for the eCQM.
9. Patients Passing Count – Displays the current number of patients meeting expectations over the total number of patients in the test deck for the eCQM.
10. Test Coverage – Displays the percentage of the measure logic that has evaluated to true for the patient test deck. This measure provides a method for determining how much of the logic has been tested.
11. Add Patient – Allows the addition of a new patient to the test deck for this eCQM.
12. Failing Patient – Displays an example of a patient that is not currently meeting expectations for the eCQM.
13. Passing Patient – Displays an example of a patient that is meeting expectations.

14. Patient Name – Displays the name given to the patient.
15. Patient Status – Displays PASS or FAIL to indicate if the patient is meeting expectations.
16. Expand Patient Results Button – Displays the calculation details of a patient. This display will show the expected and actual values for the patient against the eCQM.
17. Show Result – If a patient has been selected with the Expand Patient Results Button, this displays the result of the logic evaluated on the selected patient.
18. Show All Results – If a patient has been selected with the Expand Patient Results Button, this displays the result of each eCQM logic block evaluated on the selected patient.

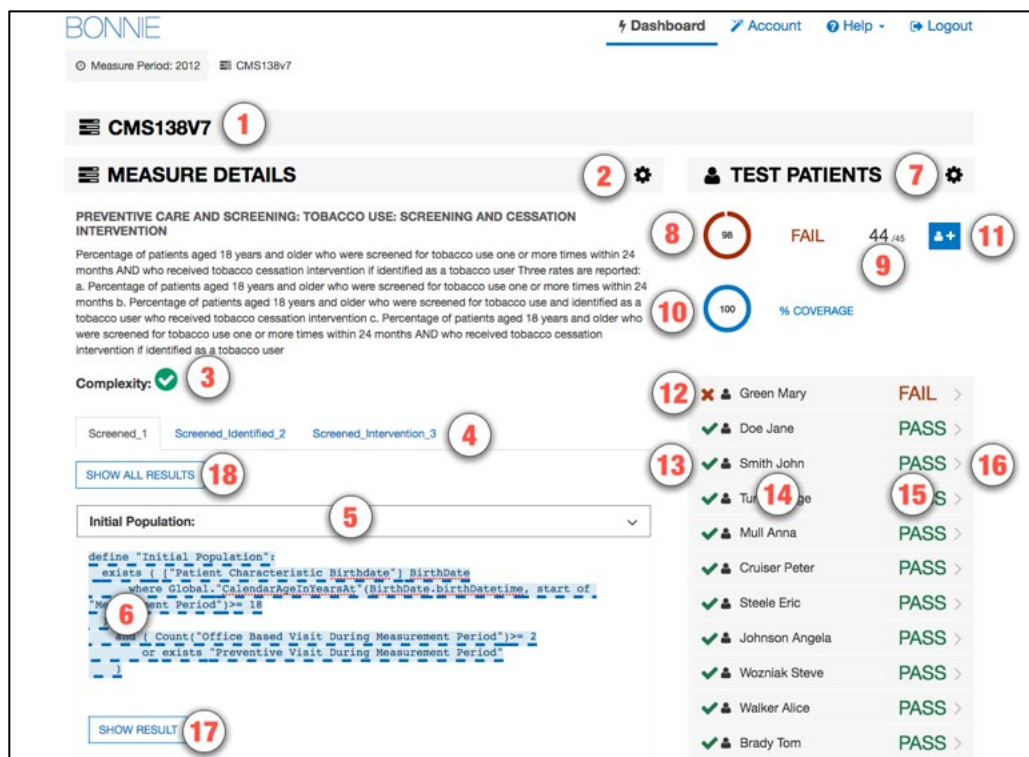


Figure 11. Measure View

The following UI elements (indicated by their item numbers in Figure 12) appear at the bottom of the Measure View Page:

1. Terminology – Displays the value sets from the eCQM and referenced libraries. Clicking on any one of these value sets will show the list of codes associated with that value set.
2. Overlapping Value Sets – Displays value sets in the eCQM and referenced libraries that have shared codes.

1 TERMINOLOGY		
Antipsychotic Medications	2.16.840.1.113883.3.464.1003.196.12.1255	Draft
Bipolar Disorder	2.16.840.1.113883.3.67.1.101.1.128	Draft
Birth date	Direct Reference Code	N/A
Dead	Direct Reference Code	N/A
Emergency Department Visit	2.16.840.1.113883.3.117.1.7.1.292	Draft
Encounter Inpatient	2.16.840.1.113883.3.666.5.3001	Draft
Ethnicity	2.16.840.1.114222.4.11.837	Draft
Huntington's Disease	2.16.840.1.113883.3.464.1003.105.12.1032	Draft
ICU Admission or Transfer	2.16.840.1.113883.17.4077.3.2040	Draft
Inpatient Psychiatric NDNQI	2.16.840.1.113762.1.4.1098.3	Draft
ONC Administrative Sex	2.16.840.1.113762.1.4.1	Draft
Payer	2.16.840.1.114222.4.11.3591	Draft
Race	2.16.840.1.114222.4.11.836	Draft
Schizophrenia	2.16.840.1.113883.3.464.1003.105.12.1104	Draft
Threat to themselves or others	2.16.840.1.113883.3.464.1003.195.12.1020	Draft
Tourette's Syndrome	2.16.840.1.113883.3.464.1003.105.12.1030	Draft
2 OVERLAPPING VALUE SETS		
Encounter Inpatient 2.16.840.1.113883.3.666.5.3001	Inpatient Psychiatric NDNQI 2.16.840.1.113762.1.4.1098.3	1 code
Inpatient Psychiatric NDNQI 2.16.840.1.113762.1.4.1098.3	Encounter Inpatient 2.16.840.1.113883.3.666.5.3001	1 code

Figure 12. Measure Terminology and Overlapping Value Sets

4.2 Measure Logic

The left-hand side of the Measure View contains the measure title, description, and a representation of the logic. The representation of the eCQM logic is similar to the human-readable display for the eCQM provided in the MAT measure exports.

The eCQM logic section can be used to verify that the eCQM logic was properly loaded from the HQMF without the loss of any logical conditions. The logic can also be used to evaluate the nature of the calculation of a test patient against the logic (please refer to subsection 4.4, *Calculation Results*) and to visualize the test coverage of the eCQM logic.

4.3 Creating a New Test Record

To create a new test record, begin by clicking the “Add Patient” button (item #9) in the Measure View (Figure 11). This action opens the Patient Builder (shown in Figure 17). After creating the test patient record, the application returns the user to the Measure View where the user can evaluate the results of calculating the patient against the eCQM.

4.4 Calculation Results

The user can calculate the test patient against the logic of the eCQM in the Measure View (Figure 11). High-level results are calculated automatically when the Measure View is loaded. These high-level results appear under the test patients section on the right-hand side of the Measure View. These results include the percent of patients passing (item #6), the test coverage (item #8), individual statuses for each patient (items #10–13), and the passing patient count (item #7).

In addition to these high-level results, detailed results can be displayed for an individual patient by clicking the “Expand Patient Results” button (#16). This displays a table of expected and actual results for the patient covering each population listed in the measure and shows whether that population currently meets or fails expectations.

The Expanded Results View employs the following UI elements (as indicated by their item numbers in Figure 13):

1. Failing Population – A population for which the patient fails.
2. Passing Population – A population for which the patient passes.
3. Population Column – A list of the population types.
4. Expected Value – The user-defined expected value for the population.
5. Actual Value – The calculated value for that population.
6. Edit Patient Button – Allows editing of the selected patient.
7. Clone Patient Button – Allows cloning of the selected patient.
8. Delete Patient Button – Allows deleting of the selected patient.

Population	Expected	Actual
IPP	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
DENOM	<input type="checkbox"/>	<input checked="" type="checkbox"/>
DENEX	<input type="checkbox"/>	<input type="checkbox"/>
NUMER	<input type="checkbox"/>	<input type="checkbox"/>

Figure 13. Expanded Results View

Clicking the “Expand Patient Results” button also displays the patient results calculated against each line of logic. This information is displayed in the measure logic section of the view by highlighting each logic clause. As shown in Figure 14, Figure 15, and Figure 16, a green

highlight (accompanied by a **solid underline**) indicates a passing result for the logic calculation, while a red highlight (accompanied by a **double-underline**) indicates a failing result over the applicable lines of text.

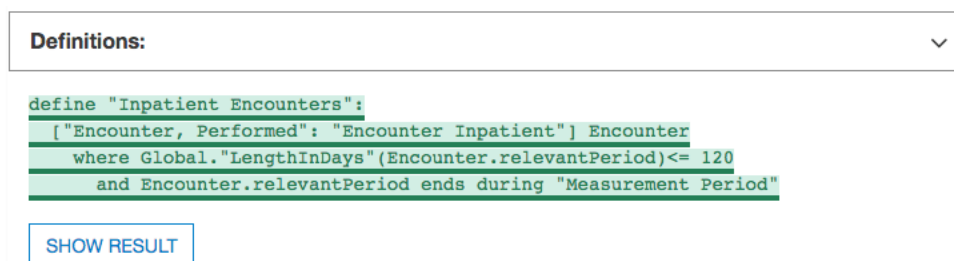


Figure 14. Logic Calculation Highlight – Passing Results

Figure 14 shows the results of a single patient calculated against the eCQM logic. The highlighting of the eCQM logic text for calculation is intended to provide a clearer understanding of why a patient is calculating to a specific population, and uses the following indications of status:

- Logic evaluated to TRUE – Green highlighting along with a solid underline
- Logic evaluated to FALSE – Red highlighting along with a double underline

The results of the calculation in Figure 14 are that the patient aligns with the logic of the initial patient population (IPP). The highlighting of the logic in Figure 14 indicates that every AND condition evaluated to true and at least one condition from each OR evaluated to true. Based on this calculation, the “Inpatient Encounters” definition evaluates to true for the patient.

The “Show Result” button at the bottom of Figure 14 can be used to expand the result, as shown in Figure 15. To collapse the results, click the “Hide Result” button.

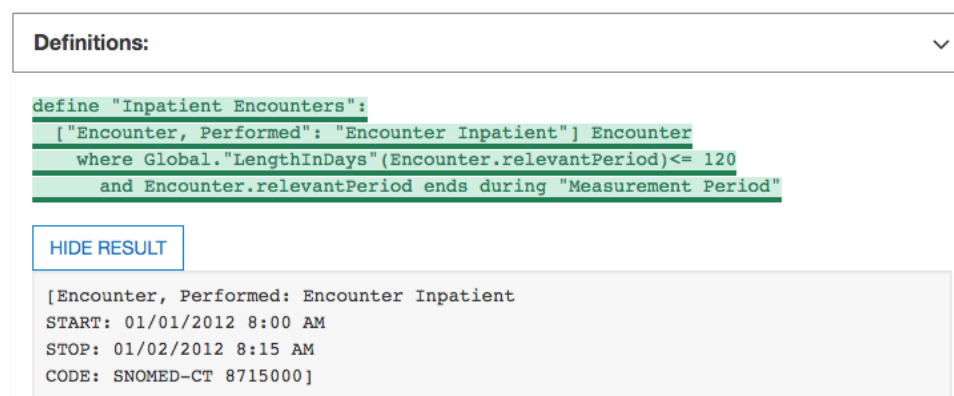


Figure 15. Logic Calculation Highlight – Expanded Result

In Figure 15, the bottom grey box displays information regarding the result of the eCQM logic definition evaluated on the selected patient. In the example shown, the “Inpatient Encounters” definition looks for encounters in a certain time range. The grey box will display all encounters

that meet that criteria; in this example, the selected patient has one encounter that meets the criteria.

In Figure 16, the population is highlighted in red, indicating that the logical statements evaluate to false. Some clauses in this figure are green, indicating they are true; however, the false female statement causes the AND logic to be false, and therefore, the results calculated in Figure 16 indicate that the patient is neither included in the IPP nor the Denominator.

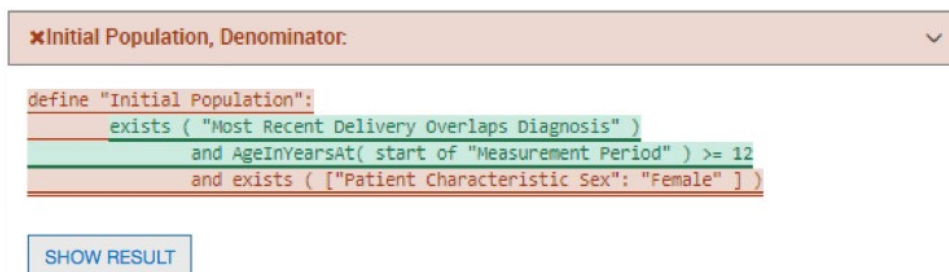


Figure 16. Logic Calculation Highlight – Failing Results

4.5 Editing a Test Record

The user can edit a test patient from the Measure View (Figure 11) by clicking the “Edit” button. The user accesses the “Edit” button for a patient (item #6 in Figure 13), accessible after expanding the patient result. Clicking the “Edit” button opens the Patient Builder (as shown in Figure 17) with the data populated for that patient. Once a patient record has been edited and saved, the application returns the user to the Measure View.

4.6 Cloning a Test Record

The user can clone a test patient from the Measure View (Figure 11) by clicking the “Clone” button (item #7 in Figure 13) to the immediate right of the “Edit” button, accessible after expanding the patient result. This action opens the Patient Builder (Figure 17) with the data populated for the patient being cloned. The difference between editing and cloning a patient is that the clone process creates a new patient based on an existing patient, while the edit process updates the data for an existing patient. Once a patient record has been cloned, edited, and saved, the application returns the user to the Measure View.

4.7 Deleting a Test Record

The user can delete a test patient from the Measure View (Figure 11) by clicking the “Delete” icon (item #8 in Figure 13) to the immediate right of the “Clone” button, accessible after expanding the patient result. Deleting a patient requires a two-step process for confirmation. Once a user deletes a patient record, the action cannot be undone. To delete a patient record, the user initially clicks the “Delete” button. A second “Delete” button is then displayed. The user must click the second “Delete” button to confirm the deletion of the patient.

4.8 Updating a Measure

The user can update a measure by clicking the “Update Measure” button, which is accessed by clicking the “Measure Actions” icon (item #2 in Figure 11). The “Update Measures” button displays the Update Measure Dialog (shown in Figure 10), which allows the user to specify a new zip package for an eCQM exported from the MAT. Once the new measure package has been entered, the user clicks the “Load” button, which updates the eCQM definition and returns the user to the Measure View with the updated measure definition.

4.9 Deleting a Measure

The user can delete an eCQM from the Measure View (Figure 11) by clicking the “Delete” icon for an eCQM. To access the “Delete” icon, the user clicks the “Measure Actions” icon (item #2 in Figure 11). A user cannot undo the deletion of an eCQM. To delete an eCQM, the user initially clicks the “Delete” icon. A second “Delete” icon is then displayed. The user must click the second “Delete” icon to confirm the deletion of the eCQM. Deleted eCQMs are no longer displayed on the Measure Dashboard.

5. Building a Patient Test Record

5.1 Overview

The Patient Builder view, as shown in Figure 17, allows the addition and editing of clinical data for a synthetic test patient record. The user accesses the Patient Builder view by clicking the “Add Patient” button (item #11) on the Measure Dashboard (Figure 5), or by clicking the “Add Patient,” “Edit,” or “Clone” buttons from the Measure View (Figure 11).

The Patient Builder View employs the following UI elements (as indicated by their item numbers in Figure 17):

1. Patient Name – Allows the entry of a first and last name for the patient record.
2. Patient Characteristics – Allows the definition of characteristics data for the patient.
3. Measure Information – Shows the description for the patient’s associated eCQM.
4. Expectations – Allows users to set the calculation expectation for each population of the eCQM.
5. Actions – Allows users to save or cancel a patient record.

Figure 17. Patient Builder View

The Patient Builder view also provides fields to either add new data or edit existing data for a patient. The patient’s first and last names can be defined using the associated edit fields. Characteristics such as gender, birthdate, race, and ethnicity can be set in the characteristics section. Quality Data Model elements can be added to the patient history section by dragging and dropping an individual element onto the patient history section.

In addition to defining the patient data, the Patient Builder view allows the user to set expectations on the patient using the “Expectations” section. Expectations represent how the user expects the patient to calculate against the eCQM. The “Expectations” section will be different based on the type of eCQM used to build the patient record. The “Expectations” section for patient-based eCQMs provides a check box for users to select the appropriate expected populations (e.g., numerator, denominator) within which the patient will be included. The “Expectations” section for episode of care-based measures allows users to select the number of episodes of care that are expected to be included in each population using a number picker. Continuous variable measures allow users to define the expected value(s) the measure is expected to calculate for the patient.

The Patient Builder’s logic section displays the logic of the eCQM against which the patient test record is constructed. As data is entered, the application continuously calculates the patient against the eCQM logic, and the results are displayed against the eCQM using the logic highlighting. Subsection 4.4 provides additional information about the descriptions of the logic highlighting technique based on calculation results shown in Figure 14 and Figure 16.


5.2 Building a Synthetic Patient

The first step in building a synthetic patient record is defining the patient name. The patient’s first and last name can be set in the “Patient Name” section (item #1 in Figure 17). The user can then add patient notes and define patient characteristics in the “Patient Characteristics” section (item #2). Characteristics for the patient include data such as birthdate, race, ethnicity, gender, payer, and living status.

After defining the patient’s name, the user defines the expectations for how the patient will behave and be calculated against the eCQM. The user sets the expectations for the patient in the “Expected Value” section (item #4). If the user expects the patient to align with the initial patient population logic, then the user should set that expectation in item #4. Expectations are used throughout the Bonnie application to determine if a patient is passing or failing against the eCQM.

If the expected results for the patient align with the actual results from calculating the patient against the eCQM logic, the patient passes. If the expected results do not match the actual calculated results, the patient fails. This system provides the capability to build patients and set expectations for those patients based on the intent of the eCQM. Once the eCQM is calculated, the passing state indicates that the eCQM logic aligns with the intent, while the failing state indicates that either the eCQM logic does not align with the intent or the synthetic patient was constructed improperly.

Note that in Figure 18, the depicted measure has a discrete expected value. Accordingly, the NUMER, DENOM, and DENEX options are shown. For continuous variable measures that rely on an aggregate calculation (e.g., a median time between ED arrival and ED departure), the OBSERV, MSRPOPL, and (optionally) MSRPOPLEX options are shown, as in Figure 18.

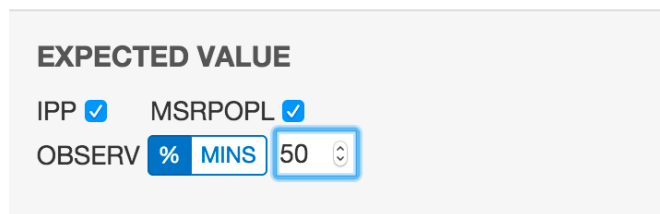

EXPECTED VALUE

IPP ☒ MSRPOPL ☒

OBSERV % MINS 50

Figure 18. Continuous Variable Measures Expected Populations

Although discrete measures offer only a checkbox for including patients in the NUMER population, continuous variable measures offer the ability to set a certain percentage or minute value for a given patient calculation. This value is only available if the patient is in the MSRPOPL population.

5.3 Building the Patient History

Figure 19 depicts four events in the patient's history—an encounter, a diagnosis, and two medications—based on data elements from the eCQM. When a QDM element is added to the patient history, it becomes an event in the patient's history, which warrants a duration and associated fields. When an event is first created, it is given default start and end date/times and is associated with a code from each value set associated with the source QDM element. These defaults as well as additional data can be edited by expanding the event. Figure 19 shows an example of an expanded event from the patient history that can be edited. By clicking the expand/collapse details icon (item #3), the user can edit the details of the element.

The following UI elements are shown in Figure 19:

1. Elements Section – Contains QDM elements that the user can add to the patient history.
2. QDM Element – Shows a condensed summary of a QDM element.
3. Expand/Collapse Details – Allows hiding or expanding the details of an element.
4. Start Date/Time – Allows setting the start date/time for an element.
5. End Date/Time – Allows setting the end date/time for an element. The checkbox allows specifying that the end date/time is undefined (ongoing event).
6. Codes Section – Allows adding codes to the element. Note that Bonnie will automatically find an appropriate code and add it.
7. Medication or Value Section – Allows adding values to the element (i.e., laboratory result values). This element shows details for adding a medication, allowing the user to specify amount prescribed, regimen, amount filled, and the times the prescription was filled.
8. Fields Section – Allows adding fields to the element (i.e., ordinality).
9. Negation Section – Allows indicating that the element is not done with a reason.
10. Delete Button – Allows deleting an element from the patient history.
11. CQL Calculation Results – Shows the results for each of the populations.

12. Show Result – Displays the result of the logic evaluated on the patient.
13. Show All Results – Displays the result for each eCQM logic block evaluated on the patient.

The screenshot displays the 'Building Patient History' interface. On the left, a vertical sidebar labeled 'ELEMENTS' (1) contains icons and labels for 'ALLERGIES INTOLERANCES', 'CHARACTERISTIC', 'ENCOUNTERS', 'INTERVENTIONS', 'LABORATORY TESTS', and 'MEDICATIONS'. The main area is titled 'PATIENT HISTORY' and shows a list of patient events. The first event is an 'ENCOUNTER' (2) dated '01/11/2012 - 01/15/2012'. Below it is an 'ALLERGY INTOLERANCE' dated '01/15/2012'. The third event is a 'Laboratory Test: Performed: LDL-c' (3). This event is expanded to show a 'RELEVANT PERIOD' (4) from '07/10/2012 8:00 AM' to '07/10/2012 8:15 AM' (5). Below the period, there is a 'CODES' section (6) with 'LOINC: 13457-7' and a 'Code Set' dropdown. The 'RESULT' section (7) shows a 'Coded' dropdown. The 'FIELD VALUE' section (8) shows a 'Coded' dropdown. The 'REFERENCES' section (9) shows a 'NOT PERFORMED' checkbox. At the bottom right of the event card is a red minus icon (10). To the right of the patient history, the 'CMS105V7' (11) section displays 'CQL Calculation Results: IPP 1, DENOM 1, DENEX 0, DENEXCEP 1, NUMER 0'. Below this are buttons for 'SHOW ALL RESULTS' (13) and 'SHOW RESULT' (12). The 'SHOW ALL RESULTS' button is highlighted with a green border. Below the 'SHOW ALL RESULTS' button is a dropdown for 'INITIAL POPULATION' with a green border. Below the 'INITIAL POPULATION' dropdown is a code block:

```
define "Initial Population":
  TJC."Encounter with Principal Diagnosis and Age"
```

. Below the code block is a button for 'SHOW RESULT' (12). Below the 'SHOW RESULT' button is a dropdown for 'DENOMINATOR' with a green border. Below the 'DENOMINATOR' dropdown is a code block:

```
define "Denominator":
  TJC."Ischemic Stroke Encounter"
```

. Below the code block is a button for 'SHOW RESULT'. Below the 'SHOW RESULT' button is a dropdown for 'DENOMINATOR EXCLUSION' with a red border. Below the 'DENOMINATOR EXCLUSION' dropdown is a code block:

```
define "Denominator Exclusions":
  TJC."Ischemic Stroke Encounters with Discharge Status"
  union TJC."Comfort Measures during Hospitalization"
```

. Below the code block is a button for 'SHOW RESULT'. Below the 'SHOW RESULT' button is a dropdown for 'NUMERATOR' with a red border.

Figure 19. Building Patient History, Including Edit Clinical Element View

After the user defines the patient characteristics, the user builds a patient history from QDM elements extracted from the eCQM. The available elements from the eCQM are organized by category and listed in the “Elements” section (item #1). The user can click a category to expand the list of available elements. Expanding this list allows the user to click and drag an individual element onto the patient history (item #2).

Several fields in the Edit Clinical Element View can be edited for an event in the patient history using the controls shown in Figure 19. These fields include the start date/time of the event (items #4 and #5), codes (item #6), values (item #7), various fields (item #8), and negation rationale (item #9). The start and end date times can be set for an event by either typing into the text fields directly or by using the date/time pickers that are displayed when the field is selected. An undefined end time can be set for the event (used for active or ongoing events) by selecting

the “Undefined” checkbox. Selecting the “Undefined” check box clears the end date/times, indicating that the event has not ended.

Values, fields, and negation rationale also can be set for an event. Values can be set by electing the type of the value (scalar or coded) and entering a scalar value or selecting a coded value from a drop-down listing of all the value sets associated with the measure. Fields such as ordinal, severity, discharge date/time, etc. can be added using the “Fields” section. Fields are added by selecting the type of the field (scalar, coded, or time) and entering a scalar value, selecting a value set, or entering a date.

Negation rationale can be defined for the event to indicate that the event was not done for a specific reason. Negation rationale is added by clicking the “Not Performed” checkbox in the negation section and selecting a value set representing the reason the event was not done.

Finally, events can be removed from the patient history by clicking the “Delete” icon (item #10). Deletion requires a two-step process. After initially pressing the “Delete” icon, the user is prompted to depress the “Delete” icon a second time to confirm the deletion.

5.3.1 Patient History Items that Are Related to Past Items

In creating the patient timeline, some items may change related to other items that occurred in the past. For example, a note or report may relate to a specific order or referral. Users can relate these two items using the “References” section, as depicted in Figure 20.



Figure 20. References Section of the Patient History Builder

5.4 Incremental Calculation

The final section of the Patient Builder View is the logic section. This section displays a representation of the logic for the eCQM against which the patient is constructed. The logic can be a reference to help describe the details of the data that should be added in building the synthetic patient.

The logic section continuously displays the results of calculating the patient against the eCQM by means of the logic highlighting described in Figure 14 and Figure 16. Any modification made to a patient triggers a recalculation of the patient against the eCQM, which updates the results of the calculation displayed by the logic highlighting. Therefore, while the user constructs the patient record, the user can inspect the behavior of the logic relative to the patient.

When the user has completed constructing a synthetic patient, the user clicks the “Save” button in the Patient Builder View. This action adds the patient to the test deck for the measure and returns the user to the Measure View. After the user creates the first patient, additional patients can be created from scratch or cloned from existing patients to extend the coverage of the test deck against the eCQM.

6. Feedback and Support

An issue tracker and feedback email list are available to support the resolution of issues and to answer questions related to the Bonnie application. The Bonnie issue tracker is available on the ONC Jira system at: <http://jira.oncprojecttracking.org/browse/BONNIE>

When users encounter bugs in the Bonnie application, they should use the issue tracker to report bugs, ask questions, or to request new features. To add an issue, users must create a login account in the Jira system. Once an issue has been entered, the Bonnie team will review and prioritize it. Alternatively, questions can also be addressed to the Bonnie feedback list at bonnie-feedback-list@mitre.org. The Bonnie feedback list email can be accessed using the “Contact” link in the main Bonnie navigation menu at the top of every page.

In addition to feedback provided through the issue tracker, users may join the Bonnie user group for community support. This group offers the opportunity for Bonnie users to ask and respond to questions about Bonnie and receive feedback from other community members. Bonnie developers will monitor and periodically provide input as well. To join the Bonnie user group, users may click on the “User Group” link on the splash page shown in Figure 21 or the User Group option in the Help menu in the application header shown in Figure 22.



Figure 21. User Group Link on Bonnie Splash Page

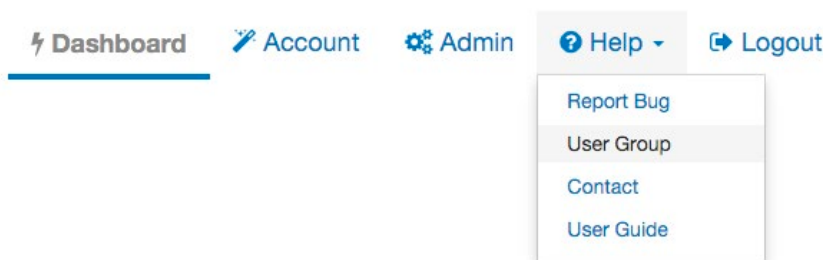


Figure 22. User Group Link in the Application Header

7. Frequently Asked Questions

Does Bonnie replace Cypress or is it an alternative to Cypress for certification?

No. Bonnie is a testing tool for measure developers to test eCQMs as they are being authored, while Cypress is the CMS quality reporting program certification tool. Bonnie cannot be used for CMS quality reporting program certification for vendors.

Bonnie provides insight into the details of how eCQM logic behaves. This makes it useful to vendors, measure developers, and others in the eCQM space who would like to better understand the behavior of eCQMs. The Bonnie testing tool uses the same measure calculation engine as the Cypress tool. Therefore, Bonnie can be used to better understand the behavior of the eCQMs when implementing an eCQM prior to certification.

Can I export patient records from Bonnie?

Users can export test patients constructed using Bonnie in a human-readable (HTML) format and in the QRDA Category 1 format.

Can I load patient records into Bonnie?

Currently, Bonnie does not support loading patient records into the tool. If you would like to calculate eCQMs using existing patients, the popHealth tool may be a better solution.

Does Bonnie automatically generate patient records?

Currently, Bonnie does not offer the capability to automatically generate patient records based on the eCQM logic. The goal of Bonnie is to allow the construction of carefully considered patient scenarios. These synthetic patient records are intended to have an expected result based on the intent of the eCQM. Automatically generating patient records would generate patients that align with the eCQM logic as written, but if there were issues in the eCQM logic, these patients might not align with the intent of the eCQM.

My patient does not match the logic of the Initial Population. Why is the patient passing?

A test in Bonnie is based on assigning an expected outcome to a synthetic patient record based on the intent of the eCQM. A patient can be constructed with the expectation that the patient calculated is to be included in the Initial Patient Population. If the patient actually calculates to the Initial Patient Population, then the test passes. If the patient does not calculate to the Initial Patient Population, then the test fails.

There must be an exact match between the patient's alignment and the expectation's needs. For instance, if the expectation is set for only the Initial Patient Population, and the patient calculates to the Initial Patient Population and the Denominator, the test fails. Similarly, a patient can be constructed with the expectation that the patient does not align with any of the eCQM populations. In this case, if the patient does not calculate into the Initial Patient Population, the test will pass. If a patient is constructed with the expectation set of not aligning with any of the eCQM populations and that patient calculates into the Initial Patient Population, that test will fail.

Where can I get help with Bonnie?

The Bonnie application has a Help menu in the header. Once you log into the application, you can use the help menu to send an email to the Bonnie feedback list to (a) ask a question directly; (b) access the Bonnie issue tracker to enter a bug, ask a question, or request a feature; and (c) access the user guide that provides step-by-step instructions on using the Bonnie tool.

What measure formats can I load into the Bonnie tool?

The Bonnie application can load CQL 1.3 measure packages from the Measure Authoring Tool using the QDM 5.4 model. Note that when loading, you will need a National Library of Medicine (NLM) Value Set Authority Center (VSAC) account (<https://uts.nlm.nih.gov/license.html>) to download the value sets associated with the measure.

Do I have to be a measure developer to use the Bonnie tool?

No. Anyone can sign up for a Bonnie account using the “register” link on the login page.

Do I need to be a Measure Authoring Tool user to use the Bonnie tool?

No. eCQMs can be loaded into the Bonnie tool either by using a Measure Authoring Tool zip file export or by loading measures released to the CMS website. Loading eCQMs using a MAT zip file export will require a MAT account to download the export zip; however, if you do not have a MAT account, you can the released versions of the eCQMs from the eCQM Library page on the CMS website once they are published.

Can Bonnie be used to calculate the results for a large number of patient records?

Bonnie is not designed to handle calculations for more than a few hundred patient records per measure. If you are interested in calculating eCQMs against a larger number of patient records, the popHealth tool may be a better solution.

Acronyms

Acronym	Definition
CMS	Centers for Medicare & Medicaid Services
CQL	Clinical Quality Language
CQM	Clinical Quality Measure
eCQM	Electronic Clinical Quality Measure
ED	Emergency Department
EH	Eligible Hospital
EP	Eligible Professional
HHS	Department of Health and Human Services
HQMF	Health Quality Measure Format
HTML	Hypertext Markup Language
IPP	Initial Payment Population
MAT	Measure Authoring Tool
NLM	National Library of Medicine
ONC	Office of National Coordinator for Health Information Technology
QDM	Quality Data Model
QRDA	Quality Reporting Document Architecture
UI	User Interface
UMLS	Unified Medical Language System
VSAC	Value Set Authority Center
XML	Extensible Markup Language